

## kVue<sup>™</sup> Stradivarius<sup>™</sup> SBRT Insert RT-4558KV1

The Stradivarius SBRT system is a complete solution for hypofractionated non-coplanar treatment. The hour glass design allows for unobstructed treatment paths. Bridges are available in two sizes and are adjustable in height. The abdominal compression paddle is also fully adjustable at the touch of a button to reduce thoracic movement.

- Compatible with the kVue System
- Indexing every 70 mm
- Laser alignment scales
- Central cut-out for unobstructed treatment paths including Non-Coplanar fields
- Optimized for use with a large bore CT scanner









Single motion lever used to secure and release the bridge.



38 cm & 48 cm Bridge RT-4558KV2 & RT-4558KV3

Bridges come in two sizes, 38 cm and 48 cm, and are completely adjustable in height. They have an easy single motion lever. They have both a horizontal and vertical scale to align easily with lasers. Their construction is completely MRI compatible.



RT-4558KV4-01



Incremented lateral placement along bridge allows for precise positioning



Release button used to raise and lower the compression paddle.

Abdominal compression paddle is fully adjustable at the touch of a button. In order to ensure compatibility with a wide variety of CT bores, the SBRT kit includes both the small and large compression paddles.

RT-4558KV2	38 cm Bridge - Adjustable Height
RT-4558KV3	48 cm Bridge - Adjustable Height
RT-4558KV4-01	Compression Paddle - Low Profile

Compression Paddle



Vac-Q-Fix<sup>™</sup> Whole Body Cushion RT-4517-230120T80F90 The SBRT package includes a 230 cm long whole body Vac-Q-Fix" cushion that is flared on one end for arm immobilization. See page 108-111 for more size options and detailed information on the full line of Vac-Q-Fix cushions.



Vac-Q-Fix<sup>™</sup> Thigh Compression Cushion

RT-4517-4079VF15



The ArmShuttle Elite is included in the SBRT package, increasing comfort and reproducibility.



RT-4489KW & RT-4489FB SofTouch™ Knee Wedge & Foot Block



Vac-Q-Fix Thigh Compression Cushions are securely held in place by non-metallic clips.





RT-4558CB01 Stradivarius™ Compression Belt

Thigh Compression Vac-Q-Fix Cushion, when combined with the SofTouch Knee Wedge & Foot Block, provides complete leg immobilization.




The Stradivarius SBRT Overlay is used on top of existing treatment and CT couch tops. It forms a base to mount SBRT bridges and compression paddles.

- Allows immobilized patient transport between CT and treatment rooms
- Optimized for use with a large bore CT scanner
- 158 KG (350 lb) patient transport load rating
- Integrated transport handles
- Indexed every 70 mm
- Laser Alignment Scales

## Stradivarius™ SBRT Overlay

Aluminum Equivalence: ~0.6 mm @ 100kVp Water Equivalence: 6 mm @ 6 MV Weight Limit: 249 kg (550 lb)

Uniformly Distributed Load

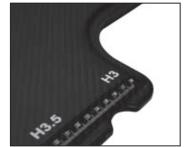

Laser Alignment Scales with indexing every 70 mm



Central Cut-Out for Unobstructed Treatments paths including Non-Coplaner Fields.